

Comprehensive selection of transducers

Relentless Innovation
for your diagnostic confidence

SAMSUNG

Curved array transducers



C2-8
Abdomen, obstetrics, gynecology, pediatric, musculoskeletal, vascular, urology



CA2-8AD
Abdomen, obstetrics, gynecology, pediatric, musculoskeletal, vascular, urology, emergency



CA4-10M
Pediatric, vascular, abdomen, obstetrics, gynecology, musculoskeletal, urology, emergency

Linear array transducers



LA3-16AD
Small parts, vascular, musculoskeletal, abdomen, obstetrics, gynecology, pediatric, emergency



LN5-12
Small parts, vascular, musculoskeletal, abdomen, obstetrics, gynecology, pediatric



EVN4-9
Obstetrics, gynecology, urology



EV2-10A
Obstetrics, gynecology, urology

Endocavity transducers

Volume transducer



VN4-8
Abdomen, obstetrics, gynecology, musculoskeletal, pediatric, vascular, urology

CW transducer



DP2B
Cardiac, vascular

Phased array transducers



PN2-4
Abdomen, cardiac, vascular, pediatric, TCD, emergency

* This product, features, options, and transducers are not commercially available in all countries.
* Due to regulatory reasons, their future availability cannot be guaranteed. Please contact your local sales representative for further details.
* This product is a medical device, please read the user manual carefully before use.

1. Optional feature which may require additional purchase.
2. Strain value for ElastoScan+™ is not applicable in Canada and the United States.
3. It is assumed that the average working hours per day is 8 hours. This is based on internal study and compared to the previous model, HM70A. The 75 hours of battery life is the combined battery life of the system and extended power pack.

SAMSUNG MEDISON CO., LTD.

© 2021 Samsung Medison All Rights Reserved.
Samsung Medison reserves the right to modify the design, packaging, specifications, and features shown herein, without prior notice or obligation.

Mobile Excellence HM70 EVO



CT-HM70 EVO V1.0-GI-eda-ZT0304-EN

CE 0123

Samsung Healthcare Cybersecurity

To address the emerging need for cybersecurity, Samsung provides a solution to support our customers by offering the tools to protect against cyberthreats that may compromise invaluable patient data and ultimately degrade the quality of care.




Mobile Excellence

The HM70 EVO ultrasound system is an impeccable example of high-performance hand-carried ultrasound system, evolved to support a diverse range of applications and patients. The HM70 EVO has exceptionally prompt workflow, durability, and high resolution imaging that can be used in a variety of clinical situations. Samsung has harmonized the key aspects of the HCU to meet today's changing needs and provide healthcare professionals the confidence they need in their work environment.



Versatility
Supports various applications and patient types

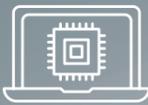

Stability
Stable performance migrating Samsung cart-based imaging engine



Intuitiveness
Intuitive interfaces for enhanced usability



Durability
Stiffness improved by 38% for various environmental usage



Mobility
Go almost a whole working day on a single charge




Imaging with exceptional quality


Exceptional image clarity and penetration is adapted from the imaging engine migrated from cart-based Samsung ultrasound systems to produce clear, uniform, high resolution 2D, 3D and color images. The HM70 EVO delivers sophisticated imaging technologies to empower ultrasound professionals with diagnostic confidence wherever needed.

Noise reduction filter to improve 2D image quality

The noise reduction filter improves edge enhancement and creates sharp 2D images for optimal diagnostic performance. In addition, ClearVision provides application-specific optimization and advanced temporal resolution in live scan mode.

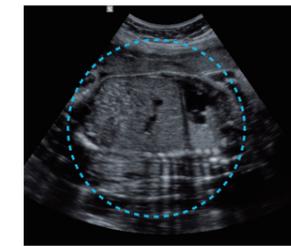


Fetal abdomen

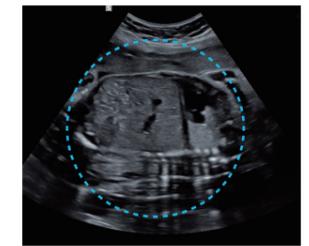


Fetal abdomen with ClearVision

Uniform imaging performance of overall image area from near-to-far

S-Harmonic™ mitigates the signal noise, enhances contrast, and provides uniform image performance of overall image area from near-to-far.

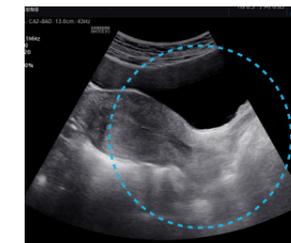

Uterus

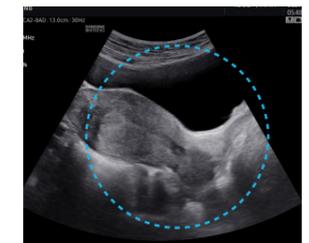

Uterus with S-Harmonic™

Directional power Doppler to examine peripheral vessels

The function uses directional power Doppler technology, enabling you to examine even the peripheral vessels. It displays information on the intensity and direction of blood flow.



Fetal circulation



Umbilical cord

3D anatomy imaging technology with detailed and realistic expression¹

RealisticVue™ displays high resolution 3D anatomy with exceptional detail and realistic depth perception. User selectable light source direction creates intricately graduated shadows for better defined anatomical structures.



1st trimester



Fetal face


Increased Operational Efficiency

The simple, intuitive interfaces enable users to immediately get to work in emergency situations. The HM70 EVO is engineered to provide the most direct interfaces and smooth actions from anywhere, avoiding delays caused by configuration. Samsung has considered even the smallest details of usability to support users with the exceptional possible performance.




Intuitively tuned interfaces

Intuitive buttons and screen layouts are optimized to express information clearly and enable complex actions to be simply performed.

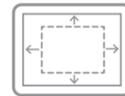

View images and cines in full screen

In the FullScreen mode, the ultrasound examination can be performed while viewing the image/cine that is fully expanded to the entire monitor.

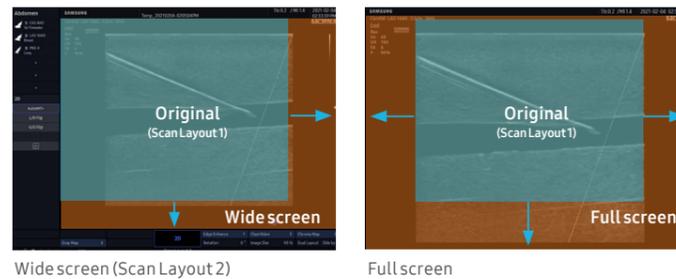

Quickly exam even in an emergency case

In case of emergency, an exam can be started with a temporary ID and patient information can be modified after the exam

Select your transducer and preset combination with one touch

With one touch, the user can select the most common transducer and preset combinations. QuickPreset increases efficiency to make a full day of scanning simple and easy.

Save image data directly to USB memory

QuickSave function allows image data to be saved directly on USB memory during the exam.

Build predefined protocols to ensure every step is followed every time¹

EzExam+™ enables users to build or to use predefined protocols. It transforms the ultrasound investigation into a streamlined process. EzExam+™ ensures the full investigation is performed, eliminating the risk of forgetting an image or loop capture, as well as measurement and transducer preset changes.

Provide valuable reference tool to help in scanning¹

As a reference tool for plane scanning, EzAssist™ has multiple clinical applications during live scanning. This allows clinicians to practice procedures using valuable reference materials such as clinical reference images and animation.

Highly Usable Design

The wave pattern of the exterior surface and the honeycomb structure of the internal design communicate strength and confidence to the healthcare professionals. The magnesium case surrounding the system and enhanced battery performance offer true usability and endurance.



Protect your system against damage from knocks and drops with the durable cover

Considering the importance of durability in a system that needs to be moved frequently, exterior wave pattern and interior honeycomb structure have been incorporated in the design to improve the rigidity of the casing by 38%*.

* Compared to the Samsung ultrasound HM70A.



Protect your system from dust, liquids, and other contaminants with KeySkin¹



1 Go almost a whole working day on a single charge^{1,3}

The battery lasts approximately 7.5 hours with the extended power pack. This is more than 100% longer than the previous model, expanding its deployment range.

2 Connect up to 3 transducers and save your time from switching over transducers¹

You can connect up to three transducers at once using the extended transducer ports on the optional cart, saving time and effort spend switching transducers in situations where every second counts.

Convenient carrying case safely protects the system with shockproof cushion¹

The carrying case allows you to safely and conveniently move your whole package including ultrasound system, transducers, and other accessories.



Available in Various Clinical Environment



Obstetrics



Gynecology



Breast




Emergency



Small Parts

The HM70 EVO supports a wide range of clinical applications, environments, and patient types—including obstetrics, gynecology, breast, emergency, and small parts. The origin, performance, and adaptability make the compact system ideal for almost any clinical situations beyond obstetrics and gynecology. The versatility enables healthcare professionals to perform effectively on targeted examinations.

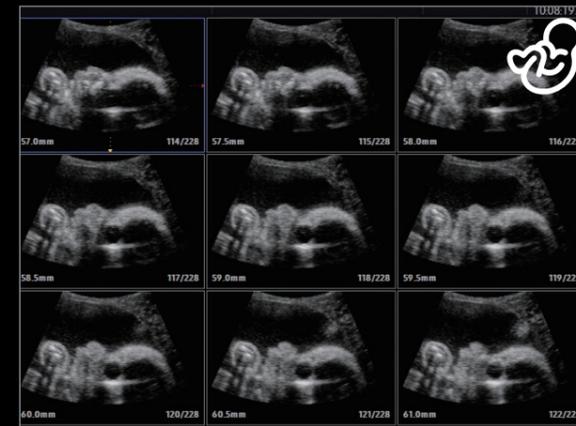


Fetal face with MSV™



Fetal heart with ClearVision



AC measurement with BiometryAssist™

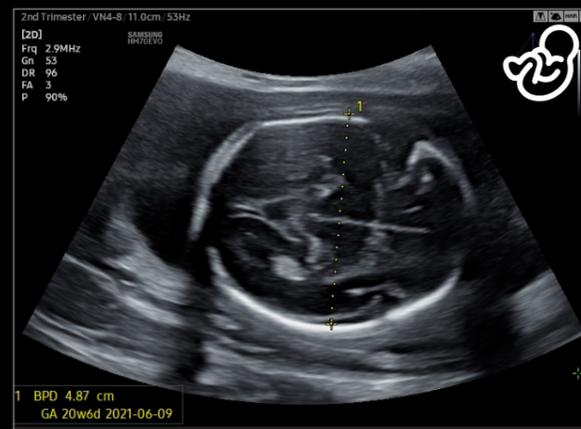


BPD measurement with BiometryAssist™




NT measurement with 2D NT¹



Fetal face with S-Harmonic™

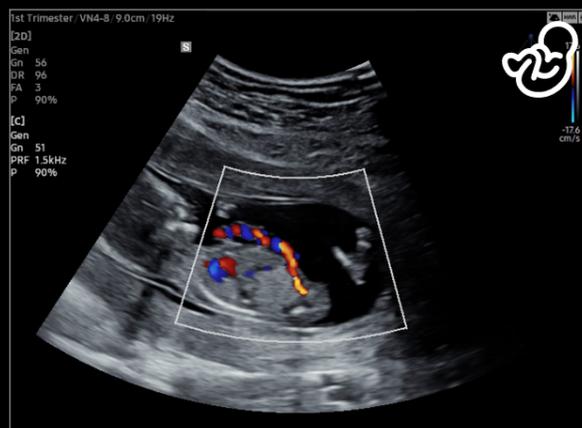

Umbilical cord

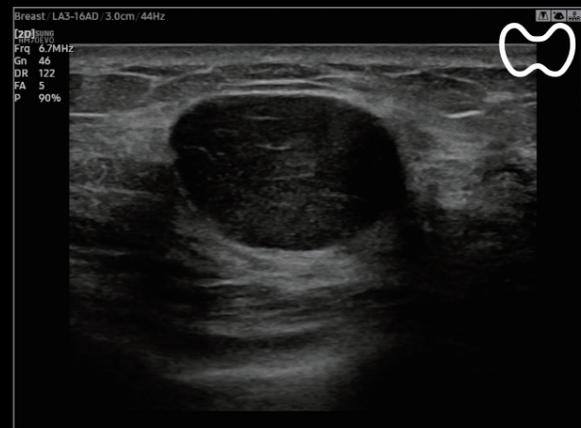

Breast mass with ClearVision



Strain ratio calculation with E-Strain™^{1,2}



Breast in 2D